


# MULTIPAD

## EFFICIENT, SIMPLE AND FLEXIBLE


The flexible MULTIPAD product family consisting of MULTIPAD Standard, MULTIPAD Plus and MULTIPAD Slim can be universally used, either alone or in combination during all examinations with head and extremity coils.

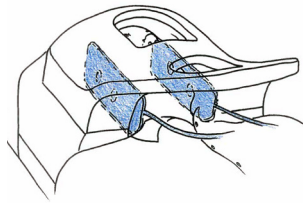



**MULTIPAD Standard, Plus or Slim**, depending on the size of the gap, are positioned on/next to headphones or directly on the patient's ears. This minimizes minute nodding and rotational movements of the head.




**MULTIPAD Slim**, either one or both pads, allows breathing movements of the patient to be reduced.

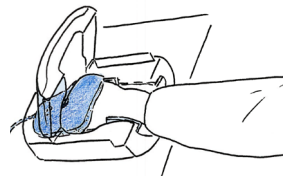

**MULTIPAD Slim**, is slid into the coil from the front and placed on the fingers. This helps reduce unwanted movements, especially due to the "Superman" position.

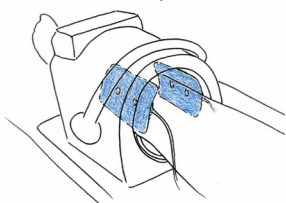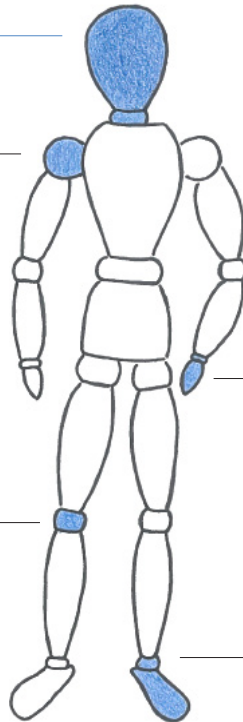

**MULTIPAD Standard, Plus or Slim** are chosen depending on the size of the gap. MULTIPAD can be positioned on the upper and lower leg laterally to the patella or used to support the stability of the lower leg.

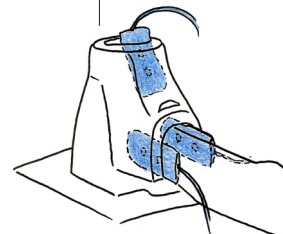

**MULTIPAD Plus**, either one or both pads are placed on the toes. This allows minor toe movement to be reduced during forefoot and upper ankle joint examinations.

Do you perform most of your examinations using flexible coils, or do you use knee and head coils for other parts of the body as well? The flexibility of the MULTIPAD product family also helps you to keep patients stable during the entire examination.

Further tips and tricks for using the MULTIPAD product family can be found on our website.