


ETCO₂ ORAL-NASAL CANNULA


LATEX FREE | SINGLE PATIENT USE | DISPOSABLE

End Tidal CO₂ Cannula simultaneously delivers oxygen while capturing CO₂. A series of strategically positioned holes deliver O₂, and CO₂ is effectively captured through the nares and scoop. This design helps lessen attenuation of the end tidal sample by oxygen dilution, particularly for mouth-breathing patients.

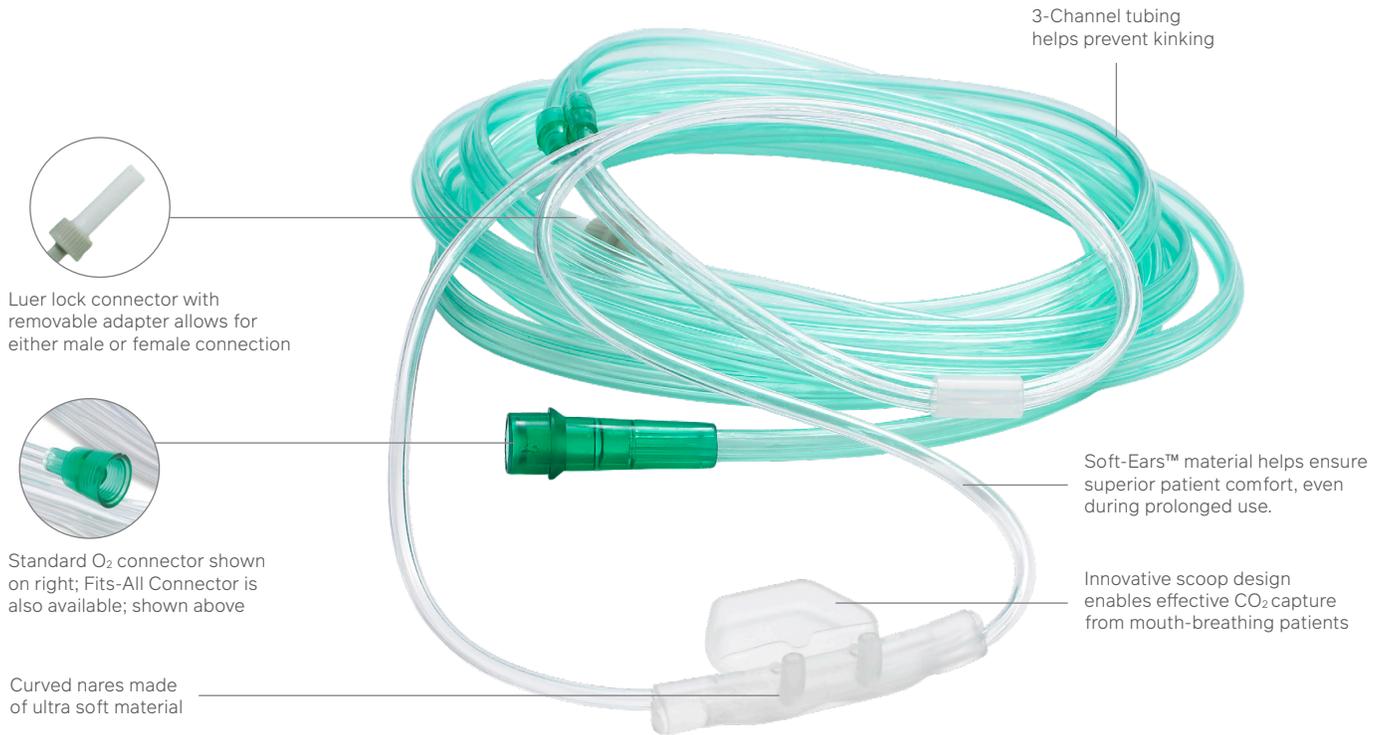

ETCO₂ ORAL-NASAL CANNULA

- Innovative design suitable for mouth-breathing patients
- Dual nasal cannula simultaneously provides oxygen to both nostrils while obtaining CO₂ sampling during spontaneous breathing
- Designed to prevent mixing of fresh oxygen with CO₂
- 3-Channel tubing helps prevent kinking
- Includes male/female adapter with luer lock connector
- Available in two tubing lengths

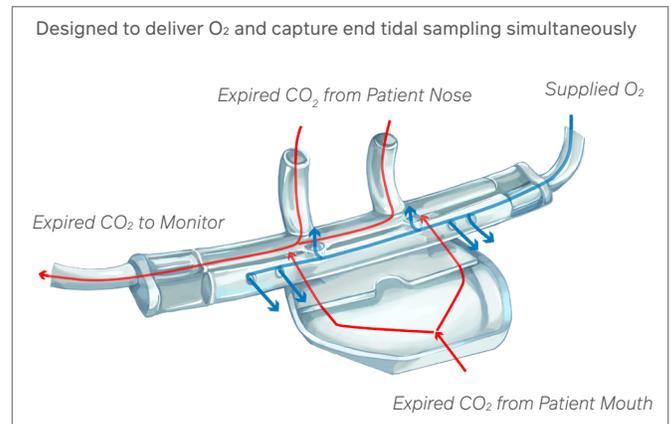

ITEM	DESCRIPTION	SIZE	LENGTH	PK
4707F	EtCO ₂ Cannula + Scoop & Fits-All Connector	Adult	7'	25

ITEM	DESCRIPTION	SIZE	LENGTH	PK
4707	EtCO ₂ Cannula + Scoop	Adult	7'	25
4714	EtCO ₂ Cannula + Scoop	Adult	14'	25

LIT-20034-C1 Rev.1